

---

# **OEM Capnography and Spirometry**

## **Products and Accessories Catalog**

Effective January 1, 2013




**CAPNOSTAT 5 CO<sub>2</sub> SENSOR**

The CAPNOSTAT 5 CO<sub>2</sub> Sensor is designed specifically for mainstream measurement of CO<sub>2</sub> using sophisticated infrared absorption spectroscopy. With mainstream, measurements are taken at the patient's airway at a sampling frequency of 100 Hz, so response is faster and there is less chance of erroneous, artifact data. The analyzer portion of the sensor cannot be contaminated with patient secretions, and there are no pumping or pneumatic components to replace. So, maintenance requirements and overall cost-of ownership are minimized. The exclusive Respironics Novamatrix family of airway adapters are easy to use and durable enough even for long-term ventilator patients. Yet, they are inexpensive and disposable. The CAPNOSTAT 5 CO<sub>2</sub> Sensor can be integrated with other components, software packages and accessories, making it easy for you to configure your patient monitoring systems with state-of-the-art CO<sub>2</sub> monitoring technology that is affordable and easy to use.

Catalog No.  
1015928

**CAPNOSTAT 5 CO<sub>2</sub> Sensor**

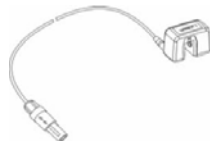

1021948

**CAPNOSTAT 5 CO<sub>2</sub> Sensor Evaluation Kit**

**CAPNOSTAT Accessories**

**SINGLE-PATIENT USE (SPU) CO<sub>2</sub> AIRWAY ADAPTERS (10 per Box)**

6063-00 SPU Adult Airway Adapters



6421-00 SPU Adult Airway Adapters w/  
mouthpiece



6312-00 SPU Neonatal Airway Adapters



**REUSABLE CO<sub>2</sub> AIRWAY ADAPTERS (Each)**

7007-01 Reusable Adult Airway Adapters



7053-01 Reusable Neonatal Airway Adapter




**CAPNO<sub>2</sub>mask O<sub>2</sub> Delivery/CO<sub>2</sub> Mainstream Monitoring Mask (10 per bag)**

9960PED-00 CAPNO<sub>2</sub>mask - Pediatric

9960STD-00 CAPNO<sub>2</sub>mask - Adult Standard

9960LGE-00 CAPNO<sub>2</sub>mask - Adult Large

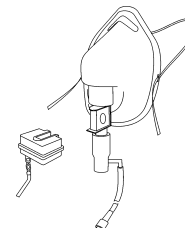

### ***LOFLO SIDESTREAM CO<sub>2</sub> MODULE***

The LoFlo CO<sub>2</sub> Module is CO<sub>2</sub> specifically designed for sidestream measurement of CO<sub>2</sub> using sophisticated infrared absorption spectroscopy. LoFlo is a sidestream sampling system with a 50 ml/minute low sampling rate that is used to measure the CO<sub>2</sub> of non-intubated and intubated patients. Because of its unique design, response is fast and there is less chance of erroneous, artifact data. The analyzer portion of the sensor cannot be contaminated with patient secretions; therefore maintenance requirements and overall cost-of-ownership are minimized. The exclusive Respirationics Novamatrix family of sidestream sampling and on-airway adapter kits are easy to use and durable enough even for long-term ventilator patients. Yet, they are inexpensive and disposable. The LoFlo CO<sub>2</sub> Module is interchangeable with the CAPNOSTAT 5 CO<sub>2</sub> Mainstream Sensor making it easy for you to configure your patient monitoring systems with state-of-the-art CO<sub>2</sub> monitoring technology that is affordable and easy to use.

Catalog No.  
1022054

Description  
***LoFlo Sidestream Module***




1040075

***LoFlo Module Evaluation Kit***

1027346

***LoFlo Sidestream Internal Engine***

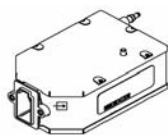

1043294

***LoFlo Internal Engine Evaluation Kit***

Catalog No.

Description

**LoFlo Accessories**

**SINGLE-PATIENT USE (SPU) CO<sub>2</sub> Cannula and Airway Adapter Kits - Without Nafion  
(10 per Box)**

3468ADU-00 Nasal CO<sub>2</sub> Sampling Cannula - Adult

3468PED-00 Nasal CO<sub>2</sub> Sampling Cannula - Pediatric

3468INF-00 Nasal CO<sub>2</sub> Sampling Cannula - Infant

3469ADU-00 Nasal CO<sub>2</sub> with O<sub>2</sub> Sampling Cannula - Adult

3469PED-00 Nasal CO<sub>2</sub> with O<sub>2</sub> Sampling Cannula - Pediatric

3469INF-00 Nasal CO<sub>2</sub> with O<sub>2</sub> Sampling Cannula - Infant

3470ADU-00 Oral-Nasal CO<sub>2</sub> Sampling Cannula - Adult

3470PED-00 Oral-Nasal CO<sub>2</sub> Sampling Cannula - Pediatric


3471ADU-00 Oral-Nasal CO<sub>2</sub> Sampling-O<sub>2</sub> Delivery Cannula - Adult

3471PED-00 Oral-Nasal CO<sub>2</sub> Sampling-O<sub>2</sub> Delivery Cannula - Pediatric

3472ADU-00 Adult Airway Adapter Set = ET >4.0mm

3472INF-00 Infant Airway Adapter Set = ET =<4.0mm

3474-00 Straight Sample Line with Luer Lok

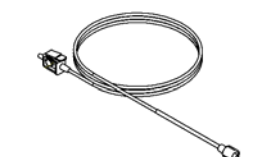

Catalog No.  
**LoFlo Accessories**

Description  
**SINGLE-PATIENT USE (SPU) CO<sub>2</sub> Cannula and Airway Adapter Kits - With Nafion  
(10 per Box)**

3468ADH-00 Nasal CO2 Sampling Cannula - Adult

3468PEH-00 Nasal CO2 Sampling Cannula - Pediatric

3468INH-00 Nasal CO2 Sampling Cannula - Infant

3469ADH-00 Nasal CO2 with O2 Sampling Cannula – Adult

3469PEH-00 Nasal CO2 with O2 Sampling Cannula - Pediatric

3469INH-00 Nasal CO2 with O2 Sampling Cannula - Infant

3470ADH-00 Oral-Nasal CO2 Sampling Cannula - Adult

3470PEH-00 Oral-Nasal CO2 Sampling Cannula - Pediatric





3471ADH-00 Oral-Nasal CO2 Sampling-O2 Delivery Cannula - Adult

3471PEH-00 Oral-Nasal CO2 Sampling-O2 Delivery Cannula - Pediatric

3473ADU-00 Adult Airway Adapter Set = ET >4.0mm

3473INF-00 Infant Airway Adapter Set = ET =<4.0mm

3475-00 Straight Sample Line with Luer Lok

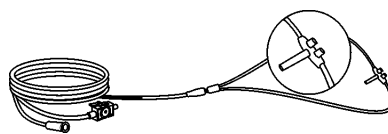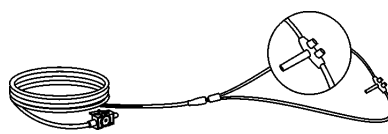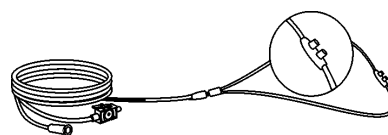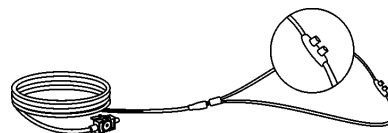

Catalog No.  
*Accessories*

Description

8751-00

Sensor Cable Holding Clips (50 per box)



1027730

LoFlo Module Mounting Bracket



1044088


LoFlo Internal Engine Cable (5 per Bag)



**FloTrak Elite Respiratory Mechanics Module**

The FloTrak Elite module measures up to twenty respiratory parameters, but is just the beginning. Combine FloTrak Elite with the CAPNOSTAT 5 Mainstream CO2 Sensor to achieve volumetric capnography and other ventilation management parameters, as well as noninvasive measurement of cardiac output, cardiac index, stroke volume and systemic vascular resistance. These parameters provide hemodynamic and ventilation management tools, guiding fluid placement and weaning protocols.

Catalog No. 1044971	Description FloTrak Elite Module
------------------------	-------------------------------------

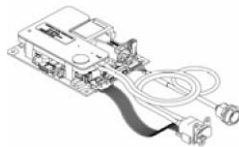

1058317	FloTrak Elite Module Evaluation Kit
---------	-------------------------------------

NOTE: The order must be entered using the FloTrak Elite P/N 1044971. An SAP sales BOM allows the evaluation kit to be ordered using the 1058317 P/N. If you order the evaluation kit P/N 1058317, the FloTrak Elite module will not be included.

Catalog No. <b>FloTrak Elite Accessories</b>	Description <b>SINGLE-PATIENT USE (SPU) Flow and Flow/CO<sub>2</sub> Sensors (10 per Box)</b>
---	--


6717-00	Pediatric/Adult Flow Sensor
---------	-----------------------------



6718-00	Neonatal Flow Sensor (box of 10)
---------	----------------------------------




6719-00	Adult combined CO <sub>2</sub> /Flow sensors
---------	--

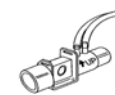

6720-00	Neonatal combined CO <sub>2</sub> /Flow sensors
---------	---



6716-00	Pediatric combined CO <sub>2</sub> /Flow sensors
---------	--

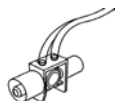







*[oem.respironics.com](http://oem.respironics.com)*

*1 (800) 243-3444. Option 3*

*(203) 697-6450*

*©2013 Respironics, Inc. All rights reserved. OEM catalog 01-January-13*